



# RELIANCE®

## ORTHODONTIC PRODUCTS

RELIANCEORTHODONTICS.COM

### Family Owned. Privately Held. Domestically Made.

ESTABLISHED 1981



MADE IN THE U.S.A.


# Reliance University



# School's in Session

First of it's kind orthodontic educational & instructional website

- ✓ Update clinical bonding techniques
- ✓ Reinforce staff training

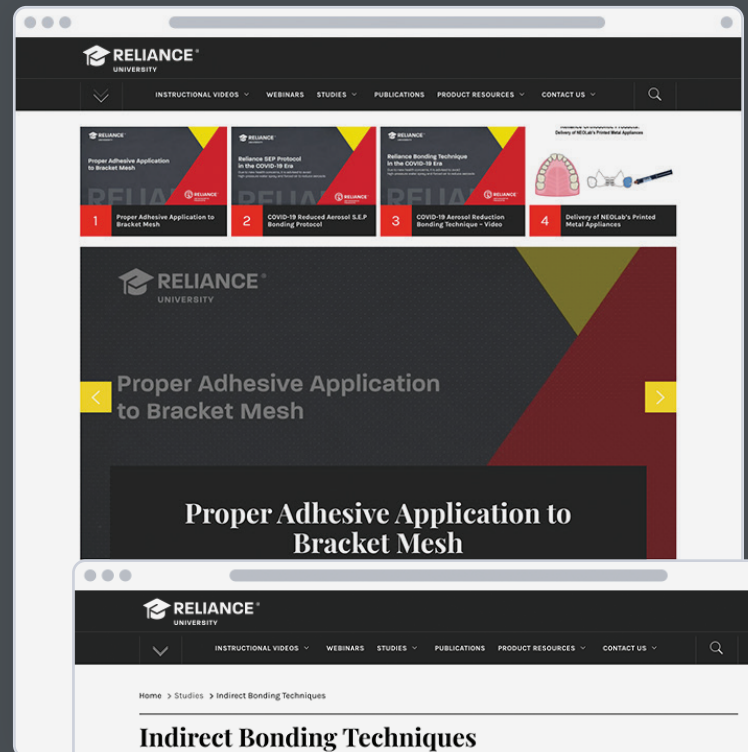

Login Now!  
RelianceU.tech

- ✓ Instructional Videos
- ✓ Literature Support
- ✓ Webinars & Livestreams
- ✓ Study Group

# Table of Contents

---

- 01 RELIANCE UNIVERSITY
- 03 HYDROPHILIC ADHESION BOOSTER
- 03 THE EVOLUTION OF ORTHODONTIC BONDING RESINS
- 04 DECAL PREVENTION
- 05 LIGHT CURE BRACKET ADHESIVES
- 09 CHEMICAL CURE ADHESIVES
- 11 SPECIALTY PROCEDURES
- 12 BOND ALIGNER
- 13 BAND/APPLIANCE CEMENTS
- 17 TOOTH PREPARATION
- 23 CLEAN UP & POLISHING
- 24 LINGUAL RETENTION
- 27 RESTORATIVE SOLUTIONS
- 28 ACRYLIC APPLIANCE BONDING
- 29 RELIAMOLD
- 30 TNT ARCHWIRES
- 31 RELASTICS
- 33 CUSTOM BASE INDIRECT BONDING
- 34 CLEAN BASE IDB ADHESIVES
- 35 ACCESSORIES
- 38 TERMS & CONDITIONS

## MADE IN THE U.S.A.

---

All Reliance adhesive products are proudly manufactured in the U.S.A.



# Hydrophilic Adhesion Booster



ASSURE® PLUS



ASSURE® PLUS UNIDOSE CAPSULES

## ASSURE® PLUS

6 ml Assure® Plus All Surface Bonding Resin	#PLUS
2 ml Assure® Plus All Surface Bonding Resin	#PLUS2ML
50 Pack Assure® Plus 0.20 ml Unidose (10 ml)	#UAP50



4 out of 5 orthodontists use Assure Plus as their self-etch or total-etch primer.

### THE MOST VERSATILE RESIN EVER DEVELOPED

- Enhances bond strength on all surfaces
- Tolerates moisture contamination
- Primes and conditions any tooth surface for bonding
- Effective on any enamel type
- Effective on any artificial substrate
- Eliminates hydrofluoric acid with porcelain crowns

### CLINICAL APPLICATIONS

- Normal or atypical enamel
- Dentin or cementum
- Eliminates additional conditioners with
  - Stainless
  - Amalgam
  - Gold
  - Composite
  - Acrylic
  - Pontics
- Eliminates hydrofluoric acid application with
  - Porcelain\*
  - Lithium Disilicate
  - Zirconia

## The Evolution of Orthodontic Bonding Resins

ADHESION BOOSTERS

ENAMEL PROTECTIVE SEALANTS



PLASTIC CONDITIONER

1980



METAL PRIMER



PORCELAIN CONDITIONER

1985



MAXIMUM CURE® SEALANT



LIGHT BOND™ SEALANT

# Decal Prevention



ORIGINAL PRO SEAL®



L.E.D. PRO SEAL®


## PRODUCT FEATURES

- Clinically proven for over 15 years
- Complete and extended enamel protection
- Full 2 year enamel protection proven and documented in numerous studies

## PRODUCT FEATURES

- Cures with any curing light
- Full 2 year enamel protection proven and documented in numerous studies

## PRO SEAL® LIGHT CURING WAVELENGTHS

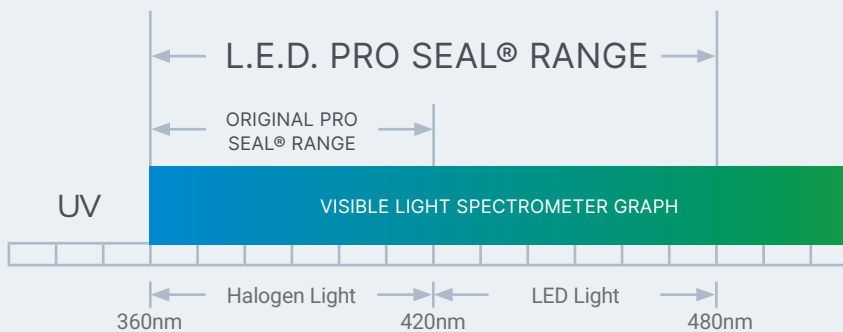

## PRO SEAL®

6 ml L.E.D. PRO SEAL®	#LEDPRO
6 ml PRO SEAL®	#PRO
UV Disclosing Light	#UVML



PLAQUE VISIBLE AT DEBOND



PLAQUE REMOVAL WITH POLISHING AT DEBOND

Pro Seal prevented plaque and tartar developing into decal by protecting the enamel integrity.



I swear by ProSeal to protect my patients' teeth from white spot lesions. The peace of mind it provides is truly invaluable!  
- Michael K. DeLuke, DDS, MDS

## INSTRUCTIONAL VIDEO

Using Pro Seal to Prevent Decalcification Available at [RelianceU.tech](http://RelianceU.tech)



## STRENGTH



ASSURE® PLUS  
Unparalleled Bond Strength

TODAY

## PROTECTION



L.E.D. PRO SEAL®  
Unparalleled Decal Prevention

# Light Cure Bracket Adhesives



PUSH SYRINGE



PUSH SYRINGE TIPS

## NEWLY DEVELOPED - SUPERIOR WORKING PROPERTIES


- Proven bond strength and consistency
- Unique viscosity easily penetrates bracket pads
- No bracket drift – stays in place until cured
- New proprietary microfiller construction
- One-piece flash removal – "ropes" off bracket periphery
- Translucent finish matches any tooth shade
- Maximum bond strength
- Easy removal and cleanup
- UV disclosing agent – fluoresces on the tooth under black light

## CLINICAL APPLICATIONS

- Bracket and tube bonding
- Fixed lingual retainer bonding
- Bonds all metal, ceramic and composite attachments
- Reliamold compatible (see Reliamolds on page 29)
- Clear aligner attachment material
  - Viscous to allow for preloading of trays

GO TO®

4g Paste Syringe	#GTP
4g Paste In Tips (10 Tips)	#GTT

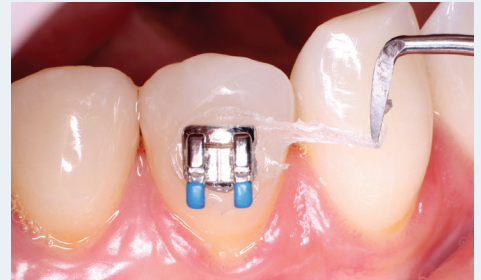

EXCESS ADHESIVE ROPES OFF IN A SINGLE PIECE FROM THE PERIPHERY OF THE BRACKET



BLACKLIGHT FLUORESCENT



CRYSTAL CLEAR FOR UNIVERSAL SHADE MATCHING WITH CLEAR ALIGNERS

RECOMMENDED PRIMER WITH GO TO® ADHESIVE

Assure® Plus (Page 3)



- \* Etchants and Conditioners appear on page 17
- \* Bonding Accessories appear on page 35
- \* Auxiliary Bonding Supplies appear on page 35



**GO TO ADHESIVE WORKING PROPERTIES ARE IDEAL FOR CLEAR ALIGNER ATTACHMENTS**



1. PUMICE



2A. ETCH



2B. ETCH



3. POST ETCH CHALKY SURFACE



4A. ASSURE PLUS



4B. ASSURE PLUS APPLICATION



4C. ASSURE PLUS AIR DRY



4D. ASSURE PLUS CURE  
(OPTIONAL BUT RECOMMENDED)



5. ADHESIVE APPLICATION



6. TRAY SEATING



7. CURE



8. FUNCTIONAL ATTACHMENT

# Light Cure Bracket Adhesives



## Multi-Purpose Composite

TKO is a versatile composite with strength and a low modulus of elasticity. TKO composite paste is effective in both laboratory and chairside applications.

ONLY FROM  RELIANCE®

### Proprietary Working Properties for Specific Applications

#### CURES STRONG AND SMOOTH WITH MODEST PLIABILITY


- Resists Opposing enamel wear
- Soft cure resin matrix for patient comfort

#### VISCOSITY PERFECTED FOR VARIED APPLICATIONS


- Controlled needle tip dispensing
- Flexes and recovers with acrylic appliances

#### PINK SHADE

- Aids in placement, removal and clean up
- Matches acrylic appliances

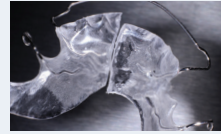


**ACRYLIC APPLIANCE REPAIRS**

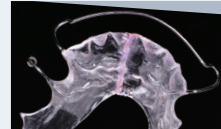


**OCCLUSAL BUILD-UPS**



**BITE OPENING RAMPS**

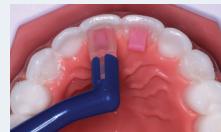

**LABORATORY BLOCK OUT APPLICATIONS**

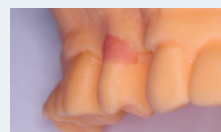

T.K.O.® (ORDER INFORMATION)

3.4g Syringe with 5 Tips

#TKOP



**INSTRUCTIONAL VIDEO**


Bonding Posterior Occlusal Buildups Available at RelianceU.tech



**PR** LIGHT BOND™  
LIGHT CURE BRACKET ADHESIVE



PUSH SYRINGE

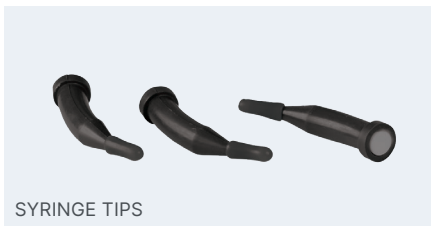

SYRINGE TIPS

**TIME TESTED PERFORMANCE**

- Unrivaled bond strength
- Highly filled paste
- Smooth tacky viscosity
- Consistent successful bonds
- Easy flash and debonding cleanup
- Designed for Doctor or staff bracket placement without drift

**CLINICAL APPLICATIONS**

- Routine bracket and tube bonding
- Intraoral bonding of metal, ceramic or composite attachments
- Indirect bonding to stone/resin models

LIGHT BOND™	SYRINGE	TIPS
Regular Viscosity Syringe*	#LBPPF	#LBPTF
Medium Viscosity Syringe	#LBMPPF	#LBMPTF
Thin Viscosity Syringe	#LBTPPF	#LBTPTF
Blue Tinted Syringe	#LBBPPF	#LBBPTF

\* Most Popular Choice



RECOMMENDED PRIMER WITH LIGHT BOND™ ADHESIVE

Assure® Plus (Page 3)

- \* Light Bond sealant appears on page 20
- Etchants and Conditioners appear on page 19
- \* Bonding Accessories appear on page 35
- \* Light Bond™ components available without Fluoride upon request
- \* All syringe applicators contain 5g of paste with Fluoride
- \* Auxiliary Bonding Supplies appear on page 29

**PR** PAD LOCK®  
LIGHT CURE FLUORESCING ADHESIVE



PUSH SYRINGE



SYRINGE TIPS

PAD LOCK®	
4.5g Pad Lock® Paste in Push Syringe with Fluoride	#PLPF
4.5g Pad Lock® Paste in Tips with Fluoride	#PLTF

**BONDS CHEMICALLY AND MECHANICALLY**

- Bonds chemically to metal bases
- Paste fluoresces under black light for flash and debonding cleanup

**CLINICAL APPLICATIONS**

- Metal bracket and tube bonds
- Indicated for fine grade mesh bases
- Non-Fluoride upon request

# Chemical Cure Adhesives

## PHASE II® MIX PASTE BRACKET ADHESIVE – CHEMICAL CURE



SCREW SYRINGE PASTE A & B



PASTE A

PASTE B

### TIME TESTED CONSISTENCY

- Microfilled paste
  - Unrivaled strength
  - Smooth, tacky viscosity
  - No bracket drift
- Self curing
  - Full polymerization
  - A&B mix system
  - 3 minute working time
- Easy removal and cleanup
- Contains Fluoride

### CLINICAL APPLICATIONS

- Routine clinical bonding
- Bracket and tube bonding

## PHASE II® DUAL CURE MIX PASTE BRACKET ADHESIVE – COMBINES CHEMICAL AND LIGHT CURING



SCREW SYRINGE PASTE A & B



PASTE A

PASTE B

### SPEED OF LIGHT CURE / SECURITY OF CHEMICAL CURE

- Initial light cure option
  - 10 second exposure tacks bracket in position
  - Allows quick and accurate positioning
- Chemical cure - time activated
  - 4 minute working time
  - Total polymerization under bases
- Contains Fluoride

### PHASE II® & PHASE II® DUAL CURE

#### PASTE COMPONENTS

Paste A 6g Syringe	#1PSAF
Paste B 6g Syringe	#1PSB
Paste A 18g Jar	#PJAF
Paste B 18g Jar	#PJB
Dual Cure Paste A 6g Syringe	#1DCSAF
Dual Cure Paste A 18g Jar	#DCJAF

#### RECOMMENDED PRIMER WITH PHASE II® ADHESIVE

Assure® Plus (Page 3)



- \* Bonding components available without Fluoride
- \* Etchants & conditioners appear on page 19
- \* Bonding accessories appear on page 35
- \* Auxillary Bonding Supplies appear on page 29

 **EXCEL®**  
CHEMICAL CURE SPECIALTY ADHESIVE



PASTE A



PASTE B

EXCEL®

Excel® Regular Paste A in 6g Screw Syringe with Fluoride	#1ERSAF
Excel® Regular Paste B in 6g Screw Syringe	#1ERSB
Excel® Regular Paste A in 12g Jar with Fluoride	#ERJAH
Excel® Regular Paste B in 12g Jar with Fluoride	#ERJB

\* Excel® Paste available w/o Fluoride on request

**TIME TESTED PERFORMANCE**

- Unsurpassed bond strength
- Thin, creamy consistency
- 3 minute working time
- Contains Fluoride

**CLINICAL APPLICATIONS**

- Indirect bonding adhesive
- Bonded acrylic appliances

 **RELY-A-BOND®**  
NO MIX BRACKET ADHESIVE



RELY-A-BOND® KIT



3.5g PUSH SYRINGE

RELY-A-BOND® COMPONENTS


3.5g Paste In Syringe	#1RPP
14g Paste In Tips	#RPT
14g Paste In Jar	#RPJ
16g Primer	#RPRF

\* All components available w/o Fluoride on request


**1 STEP / CONTACT ADHESIVE**

- 25 second working time
- Chemical cure - no light required
- Saves time and material
  - No Mixing
  - Minimizes adhesive per bond
- Highly filled paste
  - Increased bond strength
  - Tacky viscosity
  - Prevents bracket drift
- Choice of packaging
  - Push syringe
  - Syringe tips
  - Jars
- Contains Fluoride

# Specialty Procedures

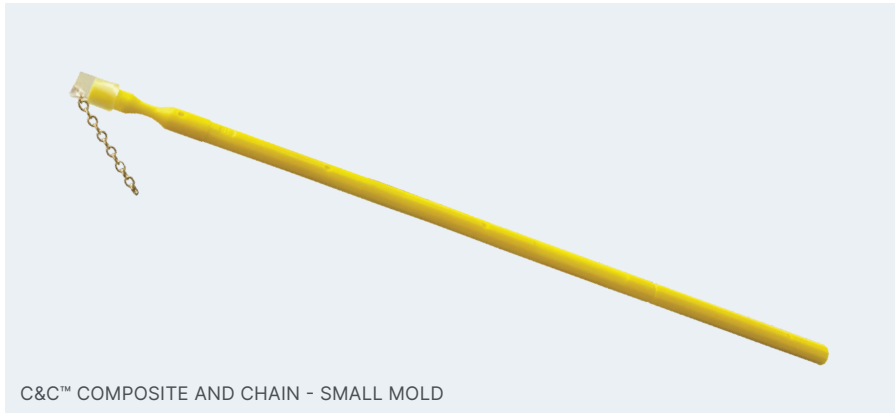

C&C™ COMPOSITE AND CHAIN - SMALL MOLD

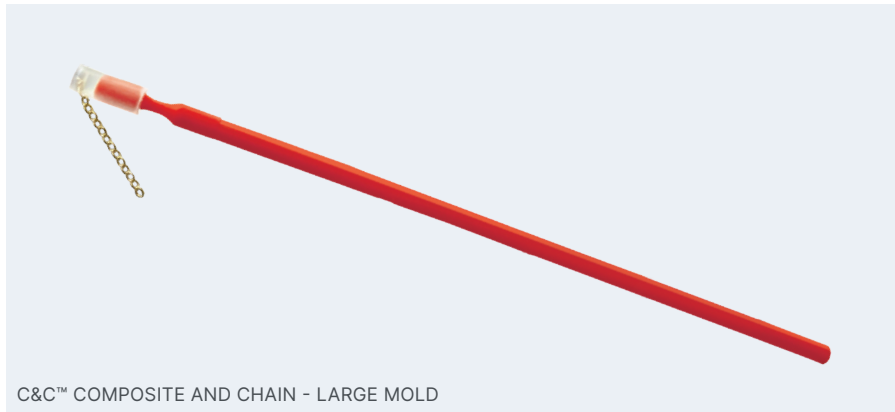

C&C™ COMPOSITE AND CHAIN - LARGE MOLD

## REVOLUTIONIZES IMPACTED EXTRUSION PROTOCOLS

- Unparalleled Strength
  - The C&C™ System does not use protruding metal attachments. Our patent pending design molds a base with composite – adapting flush to any palatal or buccal exposure.
- Painless Removal
  - Instead of painfully breaking a bracket off a mobile tooth – the C&C™ is ground down with a bur – resulting in less patient discomfort.

RECOMMENDED PRIMER WITH C&C™

Assure® Plus (Page 3)



RECOMMENDED FLOWABLE COMPOSITE WITH C&C™

FLOW TAIN™ (Page 34)




### C&C™


1 Disposable Large Attachment Mold with Chain	#DCCP
1 Disposable Small Attachment Mold with Chain	#DCCO
1.5 g FlowTain™ Flowable Composite	#FT



FILL MOLD WITH FLOWABLE COMPOSITE

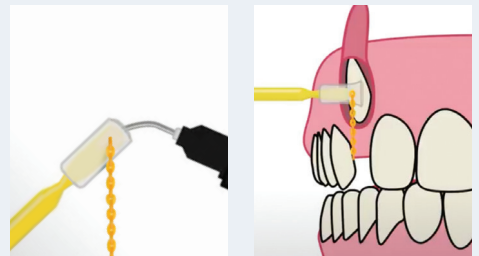

DELIVER C&C MOLD TO ENAMEL AND CURE

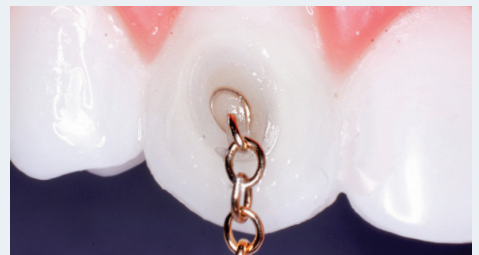

### INSTRUCTIONAL VIDEO

Utilizing Chain & Composite System Available at RelianceU.tech

# Bond Aligner



**BOND ALIGNER®**

1.4 g Bond Aligner® Syringe With 5 Tips

#BA

Button Hook Attachment

#ABUT


# Bonds Hooks and Buttons to Clear Aligners

**NO SURFACE CONDITIONING REQUIRED.**

- Will bond any button or hook to clear aligner materials - Essix, Invisalign, Invisalign Smart Track Bioplast, Tru-Tain, and more!
- Bonds to aligners without plastic preparation, primers, or conditioners.
- A clear, light-cure paste for excellent aesthetics.
- Low modulus of elasticity - Bond Aligner® flexes with thermoformed plastic trays and will not crack or separate.



## Simple Application

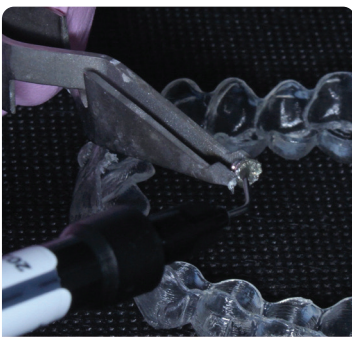

1. Apply Bond Aligner to bonding base.



2. Deliver and position.



3. Light cure.



Reliance Orthodontic Products: Over 40 Years of Serving the Orthodontic Profession.



Bonded aligner.



Bonding to the tray allows for wider clear-aligner treatment modalities

# Band/Appliance Cements



PUSH SYRINGE



ALSO AVAILABLE IN LUER LOC SYRINGE

## INDUSTRY LEADER IN BAND CEMENT


- Single Paste System
- Light Activated
- Resin base glass ionomer
  - High Fluoride release
- Natural or blue shade
- All purpose cement
  - Ultra Band-Lok® can also be used to bond large acrylic appliances and as an occlusal build-up material.
- Bonds to metal
  - Ultra Band-Lok® will bond chemically to stainless steel bands. Removal is all in one piece with cement in the appliance, making cleanup effortless.

## CLINICAL APPLICATIONS

- Band cementation
- Occlusal build ups
  - Occlusal buildup instructional video available at RelianceU.tech
- Acrylic appliances - Hyrax, MARA, RPE, etc.
  - Acrylic Appliance instructional video available at RelianceU.tech
- Lab block out material
- Non-Fluoride upon request

### ULTRA BAND-LOK®

5g Ultra Band-Lok® Syringe	#UBLPP
5g Ultra Band-Lok® Syringe Blue Shade	#UBLBPP
5g Ultra Band-Lozck® Syringe w/ Dispensing Tips	#LLUBLP
5g Ultra Band-Lok® Syringe Blue Shade w/ Dispensing Tips	#LLUBLBP
(50) Luer Loc Dispensing Tips	#TIPS50
(5) Luer Loc Dispensing Tips	#TIPS5

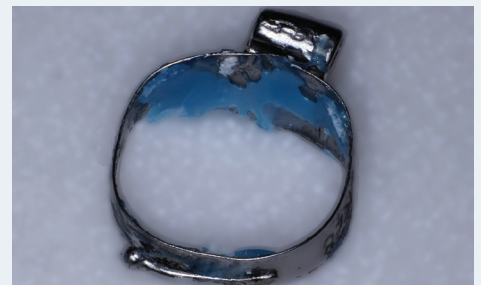

BONDS TO METAL - REMOVES WITH THE APPLIANCE AT DEBAND



OCCLUSAL BUILDUP



Bonding Posterior Occlusal Buildups  
Available at RelianceU.tech



# Revolutionary New Glass Ionomer Band Cement

- Adheres aggressively to metal
- Secure band adhesion
- Bonds in severely contaminated fields (saliva, blood, or water)
- High compressive & tensile strength
- Contains fluoride – protects enamel

ONLY FROM  **RELIANCE®**



### GLASSLOK STARTER KIT \*

- 15g Powder
- Spatulas
- Measuring Scoop
- 15cc Liquid
- Mixing Pads

### GLASS LOK® ECONOMY KIT \*

- 100g Powder
- Spatulas
- Measuring Scoop
- (3) 25g Liquid
- Mixing Pads



- Adheres to metal during removal
- Reduces cleanup time
- Resists washout under band



- Powder/Liquid mix
- Dual cure – Chemical or light activated



- Smooth consistent mix
- Natural or blue shade

**\*ALL COMPONENTS ARE  
SOLD SEPARATELY**

Order Online [www.RelianceOrthodontics.com](http://www.RelianceOrthodontics.com)

GLASS LOK®	
Glass Lok® Natural Shade Starter Kit*	#GLS
Glass Lok® Blue Shade Starter Kit*	#GLSB
Glass Lok® Natural Shade Economy Kit*	#GLE
Glass Lok® Blue Shade Economy Kit*	#GLEB
Glass Lok® Powder (15g)	#GLP15
Glass Lok® Powder - Economy Size (100g)	#GLP100
Glass Lok® Liquid Natural Shade (15g)	#GLL15
Glass Lok® Liquid Natural Shade Economy Size (25g)	#GLL25
Glass Lok® Liquid Blue Shade (15g)	#GLL15B

# Band/Appliance Cements



SELF MIX DUAL SYRINGE & SELF MIX CANULA TIPS



SCREW SYRINGES

"A" PASTE

"B" PASTE

## THE ORIGINAL, TIME-TESTED COMPOMER BAND CEMENT

- Bonds chemically to metal
  - Secure band cementation
  - Easy, clean removal
- Self mix dual syringes
  - Tips mix and dispense paste
  - Proportional paste mix
  - Eliminates spatulating
- Dual cure - light and/or chemical
  - Initial tack cure, yet full polymerization is always achieved
  - Four minute working time
- Natural or blue shade
- Contains Fluoride; available without Fluoride upon request



**#1 RECOMMENDED CEMENT FOR PRINTED METAL APPLIANCES**

### BAND LOK®

5gm Paste A / 5gm Paste B (Includes 15 Self Mix Tips) Natural Shade	#SMBL
---	-------

5gm Paste A / 5gm Paste B (Includes 15 Self Mix Tips) Blue Shade	#SMBLB
--	--------

Self Mix Tips (15)	#SMT
--------------------	------

### SCREW SYRINGES 6g PASTE


Paste A Syringe Natural Shade	#BLSA
-------------------------------	-------

Paste A Syringe Blue Shade	#BLSAB
----------------------------	--------


Paste B Syringe	#BLSB
-----------------	-------



SELF MIX SYRINGE ELIMINATES SPATULATION

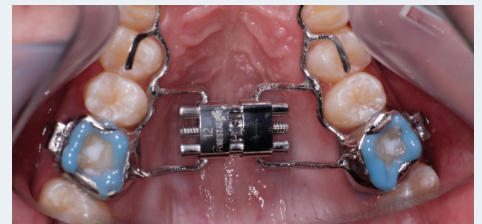


PERFECT VISCOSITY FOR PRECISE APPLICATION

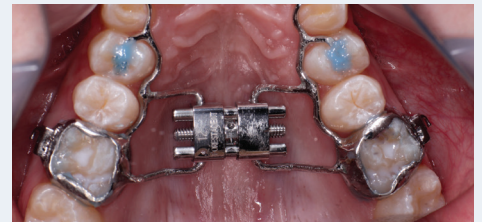

ALWAYS ACHIEVES FULL POLYMERIZATION



TACK CURED  
FLASH REMOVED  
IN ONE PIECE

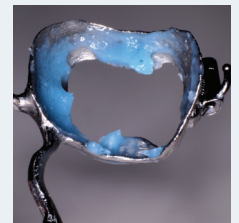

METAL BONDING  
REDUCES CLEAN UP  
AT DE-BAND



**CROWN LOK™**  
METAL CROWN LUTING CEMENT



SELF MIX DUAL SYRINGE (BLUE SHADE ONLY) & SELF MIX CANULA TIPS



SCREW SYRINGES

"A" PASTE

"B" PASTE

CROWN LOK®

5gm Paste A / 5gm Paste B (Includes 15 Self Mix Tips) Blue Shade	#SMCLB
--	--------

Self Mix Tips (15)	#SMT
--------------------	------

SCREW SYRINGES 5g PASTE

Paste A Syringe Natural Shade	#CLSA
-------------------------------	-------

Paste A Syringe Blue Shade	#CLSAB
----------------------------	--------

Paste B Syringe	#CLSB
-----------------	-------

**SECURE METAL CROWN AND BAND CEMENTATION**

- Bonds chemically to stainless steel
- Simplified removal and cleanup
- Two minute working time
- Five minute set time
- Natural or blue shade
- Contains Fluoride
- Syringe choices
  - Screw Syringe
  - Dual syringe - self mix
    - Tips mix and dispense paste
    - Proportional paste mix
    - Eliminates spatulating

# Tooth Preparation



ETCHMASTER®

STARTER KIT

Etchmaster Starter Kit	#EMHA
------------------------	-------

SMALL TIPS

Regular Tips 25 Pk	#EMTRS
Lingual Tips 25 Pk	#EMTLS
Regular Tips 50 Pk	#EMTRS50
Lingual Tips 50 Pk	#EMTLS50

LARGE TIPS

Regular Tips 25 Pk	#EMTRL
Lingual Tips 25 Pk	#EMTLL
*90 Micron 25 Pk	#EMTRL90




SURFACE MICROETCHING

**INSTRUCTIONAL VIDEO**


Bonding to Various Artificial Substrates  
Available at [RelianceU.tech](http://RelianceU.tech)

- \* Etchmaster Handpiece Starter kit includes
  - 4 Hole Handpiece Attachment
  - (5) Regular Small aluminum Oxide Tips
  - (5) Lingual Small aluminum Oxide Tips
- \* Recommended for removing composite from bracket bases and microetching stainless steel appliances. Not recommended for use on enamel.

REGULAR TIPS

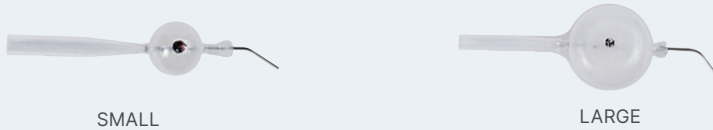

LINGUAL TIPS

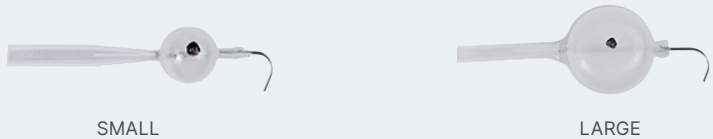

**MICROETCHING MADE CLEAN AND EFFICIENT**

- Precise surface microetching
  - Increases surface retention
  - Insures properly prepped surface
  - Maximizes bond strength on all surfaces
- 4 Hole Handpiece
  - Standard airline attachment
  - Easy chairside engagement
  - Stainless steel construction
- Preloaded cartridge tips
  - Disposable, no-mess application
  - No aluminum oxide management
  - 50 micron aluminum oxide
  - 2 tip sizes Small (6 teeth) or Large (18 teeth)

**CLINICAL APPLICATIONS**

- Lingual bonding
- Artificial substrate bonding
- Molar prep for bonds or bands



ASK®

ASK® Kit

#ASK

# One Kit All Surfaces



The new **ASK™** (All Surface Kit) from Reliance Orthodontics is all you need to bond to any intraoral surface:

- ✓ Composite
- ✓ Enamel
- ✓ Porcelain
- ✓ Zirconia
- ✓ Metal

- ✓ Reduce your bonding & inventory costs!

# Tooth Preparation



USEP®

## ONE STEP ETCHANT AND PRIMER

- Combines preparatory etching and priming into one simple step
- Effective and efficient tooth prep
- Reduces chairtime
- Single liquid application – No mixing
  - Scrub onto enamel for 5 seconds
  - Air Dry
  - Apply Assure Plus
  - Air dry
  - Deliver bracket and cure
- Effective for all enamel surfaces
  - Bracket and Tube bonding
  - Buildups
  - Increasing band cement retention
  - Reliamold delivery

INCREASE ANY S.E.P. STRENGTH WITH ASSURE PLUS


Assure® Plus (Page 3)



USEP®


6ml Bottle

#USEP

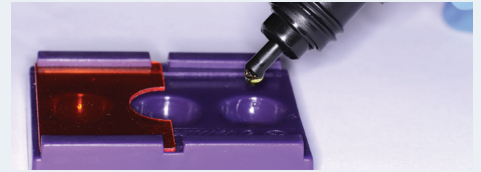

DISPENSE ULTRA S.E.P. INTO CHAIRSIDE WELL, NO MIXING REQUIRED

OR

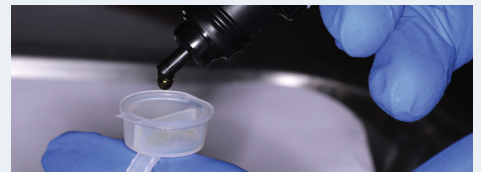

DISPENSE ULTRA S.E.P. INTO RELIANCE DUAL RING WELLS. NO MIXING REQUIRED

OR



DRASTICALLY REDUCE MATERIAL COST ON REBONDS BY SATURATING A MICROBRUSH WITH A SINGLE-TOOTH VERTICAL DISPENSE



WITH A MICRO-BRUSH, AGGRESSIVELY SCRUB ULTRA S.E.P. ONTO TOOTH FOR 5 SECONDS



AIR DRY FOR 2-3 SECONDS BEFORE PLACING ATTACHMENT/APPLIANCE

- \* Chairside well appears on page 36.
- \* Dual ring wells available soon.

 **RELIANCE S.E.P.™**  
SELF ETCHING PRIMER



RELIANCE S.E.P.™

S.E.P.™ Kit (Contains all components)	#SEPPRO
Refill Cartridges (4ml)	#SEPR
Click Syringe Dispenser	#SEPD
Mixing Well w/ Light Impervious Shield	#1SEPW

RELIANCE S.E.P.™

- Etchant and bonding primer
- One Step application
- Click-syringe dispenses pre-measured amounts
- Each dispensing click conditions 12-14 teeth
- No need to etch before applying S.E.P

 **PREP STEP™**  
DISPOSABLE PROPHY ANGLE



PREP STEP™

Prep Step™	#DPA
------------	------

\* Kits and components available w/o Fluoride on request

PREP STEP™

- 144 Individually wrapped and disposable Prep Step™ Prophy Angles are manufactured to run smooth and vibration free. Latex free, 4 web and soft cup are designed to fit a slow speed handpiece. Just throw away after use!

# Tooth Preparation

## BLUE GEL ETCHANT



BLUE GEL ETCHANT

ETCHMASTER® HANDPIECE

- Controlled and even application
- Excellent for lingual retention (page 25) or large acrylic appliances (page 28)
- 35% phosphoric acid
- Contains antimicrobial B.A.C

### BLUE GEL ETCHANT

Contains 5g syringe Gel Etchant with 10 tips #BEGS



BLUE SHADE ENSURES UNIFORM COVERAGE

## GEL ETCHANT



SYRINGE W/ TIPS 9g BOTTLE 18g JAR

- 37% phosphoric acid
- Controlled application

## LIQUID ETCHANT



LIQUID ETCHANT

- 37% phosphoric acid
- Green disclosing tint
- Regular or cherry flavor

### GEL ETCHANT

6.5g Syringe w/ 20 Tips	#EGSS
9g Bottle	#EG9
18g Jar	#EG18

### LIQUID ETCHANT

9g Bottle	#EL9
23g Bottle	#EL24
23g Cherry Flavor	#EL24C

\* Plastic Conditioner appears on page 28


## PORCELAIN CONDITIONER FOR BONDING TO CERAMIC CROWNS



PORCELAIN CONDITIONER

- Effective with Porcelain, zirconia and lithium disilicate

## PORC ETCH™ PORCELAIN ETCHING AGENT WITH BARRIER GEL

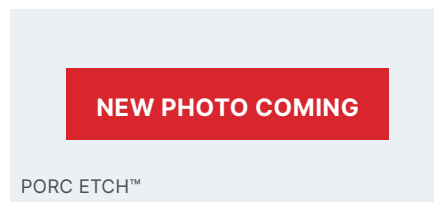

PORC ETCH™

- Preps porcelain crowns for bonding

### PORCELAIN CONDITIONER®

8cc Porcelain Conditioner	#PO
---------------------------	-----

### PORC ETCH™

7g Porc-Etch™ w/ Barrier Gel	#POE
------------------------------	------

\* Ceramic bonding instructional video available at [RelianceU.tech](http://RelianceU.tech)

 **ASSURE®**  
THE ORIGINAL UNIVERSAL BONDING RESIN

ASSURE®

6ml Assure® Universal Bonding Resin	#AS
2ml Assure® Universal Bonding Resin	#AS2ML



ASSURE®

- Bonding primer and enhancer
- Typical or atypical enamel
- Metal and composite restorations

 **LIGHT BOND™ SEALANT**  
LIGHT CURE RESIN

LIGHT BOND™ SEALANT

	REGULAR	FILLED
3g	#LBS3F	#LBFSF
7g	#LBS7F	#LBFS7F



LIGHT BOND™ SEALANT

- Bonding primer and sealant
- Available with or without quartz filler

# Clean Up & Polishing



## STEP 1 – BULK COMPOSITE REDUCTION



- 18 fluted carbide
- Quick removal of bulk composite
- Gentle on enamel
- Use with or without H2O

Carbide Bur 118S #RSB



Carbide Bur 118L #RSBL



Carbide Bur 218 #RSB2



Carbide Lingual Bur 815 #RSB815



## STEP 2 – RESIDUAL COMPOSITE REMOVAL



- Polishes tooth surface
- Black point discloses when composite is fully removed
- Removes residual composite
- Does not create heat

Straight Hand Piece (Pack of 12) #RSPHP



Latch (Pack of 12) #RSPL



Friction Grip (Pack of 12) #RSPFG



## STEP 3 – FINISHING & POLISHING



- Pre-packaged—sterile cups for individual patient use. Pre-mixed—just open and apply with prophy cup.
- No added chemicals, flavoring, fluoride or oils.
- Enamel polisher—1st & Final® can be used as final polishing step with the Reliance Renew™ Finishing System.
- Restore Diamond Polishing Paste for restoration of glaze on artificial surfaces.

1st & Final® (200 sterile cups) #FF



2g Syringe Restore™ #RESPP2





# Lingual Retention



L.C.R.™

L.C.R.™	
5g L.C.R.™ Paste in Push Syringe	#LCR
1.5g L.C.R.™ Paste in Push Syringe w/ Tips	#LCR15
5g L.C.R.™ Paste in Tips (11)	##LCRT5

- Micro+ nano fill for flowable consistency with durability of highly filled composite
- Effortlessly slumps over wire
- Fluoresces w/blacklight
- Cures to a glassy finish
- Patient comfort - reduced wear
- Developed for BLR applications



GOTO®



FLOWTAIN™

GOTO®	
4g Paste Syringe	#GTP
4g Paste In Tips	#GTT
FLOWTAIN™	
1.5 g Flowtain™ With Tips (10)	#FT
19 Gauge Tips (10)	#19TIPS10

## BLR APPLICATION


- Good working viscosity
- Increased strength
- \* Full description on page 5

## BLR APPLICATION

- Thixotropic viscosity
- Cures glassy smooth
- Slumps over wire
- Fluoresces w/blacklight
- \* Full description on page 33



ASSURE® PLUS

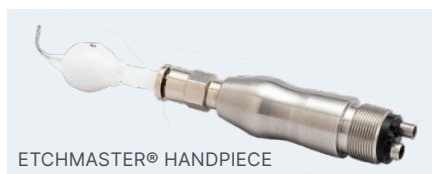

ETCHMASTER® HANDPIECE

ASSURE® PLUS	
6 ml Assure® Plus All Surface Bonding Resin	#PLUS
2 ml Assure® Plus All Surface Bonding Resin	#PLUS2ML
50 Pack Assure® Plus 0.20 ml Unidose (10 ml)	#UAP50
100 Pack Assure® Plus 0.20 ml Unidose (20 ml)	#UAP100

## BLR APPLICATION

- Enhances bond strength
- Overcomes mild contamination
- \* Full description on page 3

## BLR APPLICATION

- Mechanically prepares tooth surface
- \* Full description on page 17

ETCHMASTER® HANDPIECE	
Etchmaster Starter Kit	#EMHA

# Lingual Retention

## **ORTHO FLEXTech®** **FLEXIBLE LINGUAL RETAINER** STAINLESS STEEL OR WHITE GOLD



PATENTED

### A NEW GENERATION IN FIXED RETENTION

- Secure and stable long term passive retention
- Patented interlocking chain design
- Adapts to lingual anatomy
  - No bending or burnishing
- No models or impressions
- Flexes interproximally
  - Resists breakage and bond failures
- No reciprocal forces – allows for minor biological movement – will not take an active bend
- Stainless steel offers increased tensile and torsional strength
- Different material choices
  - White Gold (original)
  - Nickel Free Gold
  - Nickel Free Gold – etched
  - Stainless Steel – increased strength
- Low profile – patient comfort

### CLINICAL APPLICATIONS

- 3 to 3 retention
- Periodontal splinting

#### ORTH FLEXTech®

30" SS Spool	#SSOFT
30" SS Spool Etched	#SSOFTE
30" Gold Spool (white gold)	#OFT2
30" Nickel-Free Gold	#NFOFT2
30" Nickel-Free Gold-Etched	#NFOFT



PATENTED



**INSTRUCTIONAL VIDEO**

Ortho Flex Tech  
Available at [RelianceU.tech](http://RelianceU.tech)

## **EXTEND® LTR** RETENTION WITH TWO BOND POINTS



EXTEND® LTR

#### EXTEND® LTR

##### INDIVIDUAL SIZES - 5 PER PACKAGE

LTR 18mm	#EX18
LTR 20mm	#EX20
LTR 22mm	#EX22
LTR 24mm	#EX24
LTR 26mm	#EX26
Measure Device (1)	#EXMD
Extend® LTR Starter Kit	#EXK

### 3 TO 3 LINGUAL RETENTION

- Titanium molybdenum wire
- Low profile 027 TMA wire
- Can be bonded to all 6 anteriors or to cuspids only
- Flattened and serrated cuspid pads
- Preformed with cuspid offsets
- Hygiene friendly
- Nickel free

\* Extend® LTR kit includes

- (5) 18mm wire
- (5) 20mm wire
- (5) 22mm wire
- (5) 24mm wire
- (5) 26mm wire
- (25) Arch Measuring Devices

**RETAINIUM® XL**  
OPTIMUM STRENGTH BETA TITANIUM RETAINER WIRE

RETAINIUM® XL

Retainium XL – (10) 6" Wires

#XLRET



**FIXED RETENTION WITH TITANIUM MOLYBDENUM**

- Flat .027" x .016" ribbon arch
- Low profile – comfortable and functional
- Increased strength
- MRI compatible
- Easy adaptation
- Nickel free titanium construction

**CLINICAL APPLICATIONS**

- Fixed lingual retention
- Diastema maintenance
- Space maintenance
- Splinting

**RETAINIUM®**  
ORIGINAL TM ALLOY  
RETAINER WIRE



**LINGUAL TITANIUM MOLYBDENUM RETENTION**

- Low profile
- .027" x .011" ribbon arch
- Adaptable
- Flexes Interproximally
- Nickel free
- MRI compatible

**CLINICAL APPLICATIONS**

- Fixed lingual retention

**BOND-A-BRAID®**  
LINGUAL RETAINER WIRE



**CLOSE AND ACCURATE ANATOMICAL ADAPTATION**

- The dead soft wire easily adapts and is excellent for semi-permanent orthodontic splinting of the upper and lower incisors, diastema maintenance and other tooth splinting
  - Lingual retainer
  - .027" x .011" ribbon arch wire
  - 8-strand braided wire
  - Dead soft – adapts easily

RETAINIUM®

Retainium – (10) 6" Wires

#RET

BOND-A-BRAID®

Bond•a•Braid® 10 (6") Wires

#BAB2



**INSTRUCTIONAL VIDEO**

Maxillary Retention Using Retainium  
Available at RelianceU.tech

# Restorative Solutions

## PERFECT-A-SMILE® PONTIC PAINT IMMEDIATE RESTORATIONS AND SPACE MAINTENANCE



### QUICK AND ECONOMICAL AESTHETIC SOLUTION

- Light cure paste in 10 shades
- Applies smoothly to clear aligner sockets
- Translucent to match enamel
- High elasticity modulus
- Will not crack or splinter
- Bonds to any clear aligner
- Quick and economical
- Same day fix
- No doctor chairtime
- No lab fees or delays
- Complete in-office system
- Easy temporary bridge

### CLINICAL APPLICATIONS

- Eliminates flippers
- Temporary restorations
  - Congenital, trauma
  - Pre-pontic, peg laterals
  - Grafting, implants
- Space maintenance
- Aesthetic restorations
  - Single tooth
  - Partial arch

### PERFECT-A-SMILE® PONTIC PAINT

Most Popular Shades A2, B1, D3 and C4

Applicators Included With Each Jar

	4g Jar PP/Bleach	#PPBL
	4g Jar PP/A1 Shade	#PPA1
	4g Jar PP/A2 Shade	#PPA2
	4g Jar PP/A3.5 Shade	#PPA35
Vita Shade Guide References	4g Jar PP/B1 Shade	#PPB1
	4g Jar PP/C1 Shade	#PPC1
	4g Jar PP/C4 Shade	#PPC4
	4g Jar PP/D3 Shade	#PPD3
	4g Jar PP/Light Shade	#PASL
	4g Jar PP/Medium Shade	#PASM

Super Fine Micro Brushes (100) #MBSF



## PERFECT-A-SMILE® PONTIC PUTTY LIGHT CURE FORMABLE PASTE



### FINGERTIP PONTIC FABRICATION

- Light cure composite paste
- Putty-like consistency
- Facilitates tooth shape formation
- Eliminates pontics, wax & PVS
- Bonds tenaciously to stone/resin models
- Provides proper clear aligner anatomy


### CLINICAL APPLICATIONS

- Pontic formation by hand
- Bridges/fills gaps in model
- Facilitates proper aligner formation

### PERFECT-A-SMILE® PONTIC PUTTY

5g Paste In Syringe


#LBOPS

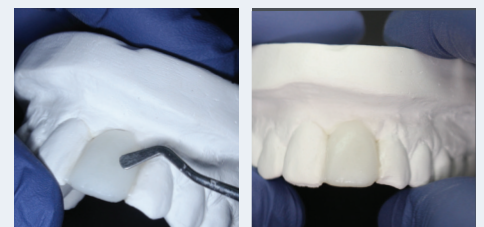

### INSTRUCTIONAL VIDEO


Creating Pontic Teeth with Perfect A Smile  
Available at [RelianceU.tech](http://RelianceU.tech)

# Acrylic Appliance Bonding

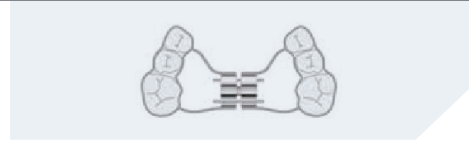

## ACRYLIC APPLIANCE PRE-DELIVERY PREP

 **BLUE GEL ETCH**  
SYRINGE APPLICATOR WITH TIPS

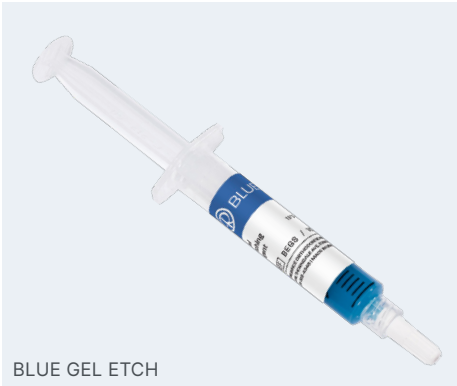

BLUE GEL ETCH

- Controlled application around buccolingual circumference
- Order #BEGS
- \* Full description on page 18

 **ASSURE® PLUS**  
BONDING PRIMER AND ENHANCER



ASSURE® PLUS

- The only primer that bonds to deciduous enamel
- Order #PLUS
- \* Full description on page 3

 **PLASTIC CONDITIONER**  
ACRYLIC PREP



PLASTIC CONDITIONER

- Initiates bond between appliance and adhesive
- 10g bottle
- Order #PL

## ACRYLIC APPLIANCE CEMENTS

 **BAND LOK® DUAL CURE**  
SELF MIX SYRINGE



BAND LOK® DUAL CURE

- Order #SMBL - Natural
- Order #SMBLB - Blue
- \* Full description on page 15

 **EXCEL®**  
CHEMICAL CURE



EXCEL®

- Order #1ERSAF & 1ERSB
- \* Full description on page 8

# Reliamold

MOLDS MAY BE DELIVERED WITH BONDER HANDLE (#BHA) OR ANY HEMOSTAT.




## Mold Composite Auxiliaries Directly on Teeth

### EASY, EFFECTIVE & ECONOMICAL


- Molds are translucent and flexible to facilitate composite loading, curing, and disengagement.
- A mold is filled with a light cure composite (Go To, TKO, or Flowtain recommended) then delivered to the tooth and cured.
- The desired auxiliary is delivered, formed, and functional immediately.
- Value priced to be single-use disposable.



EASY-TO-USE  
PRECISION MOLDS

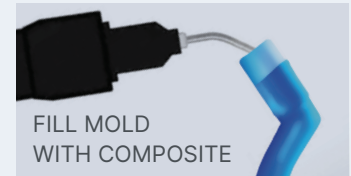

FILL MOLD  
WITH COMPOSITE

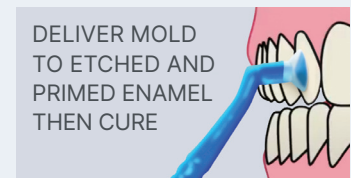

DELIVER MOLD  
TO ETCHED AND  
PRIMED ENAMEL  
THEN CURE

### DIFFERENT MOLDS FOR DIFFERENT APPLICATIONS

ONLY FROM  RELIANCE®



#### HABIT CONTROL

Inhibits thumbsucking and tongue thrust tendencies.

Habit control #RMH



#### CUSPID RAMP

Bite opening ramp applied to the cuspid—5.0 mm.


Cuspid ramp #RMC



#### BUTTON

Low profile and smooth. Facilitates elastomeric or ligature engagement.


Button #RMBU

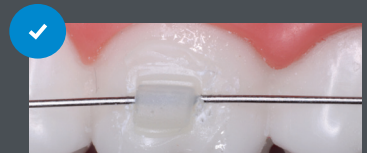

#### MOLAR TUBE

Molds a temporary tube. Accepts wire up to 020.

Molar Tube #RMT

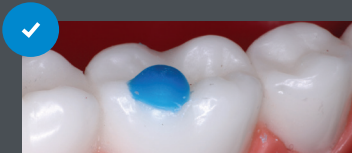

#### MOLAR PAD

Open the bite with posterior occlusal buildups. TKO or Ultra Band Lok recommended.


Molar Pad #RMMP



#### BITE RAMP

Bite opening turbo formed on the lingual of upper incisors. 2 sizes—3.0 mm or 5.0 mm.

Bite ramp 3mm #RMBR3  
Bite ramp 5mm #RMBR5

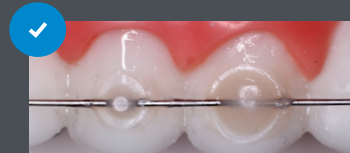

#### WIRE HOLDER

Bonds lingual retainers or other wires consistently with a smooth and low profile finish. 2 sizes—3.0 mm or 4.5 mm.

Wire bonder large #RMWBL  
Wire bonder small #RMWBS

### MOLD KIT

(2) Bonder Handles  
(108) Molds; 12 Per Type

Mold Kit #RMK

### INDIVIDUAL MOLDS

(25) Molds Per Package

Bonder Handle Sold Separately  
#BHA


# TNT Archwires



**TNT®**  
TRIED N' TRUE

## ARCHWIRES

STAINLESS STEEL / NICKEL TITANIUM  
TITANIUM MOLYBDENUM

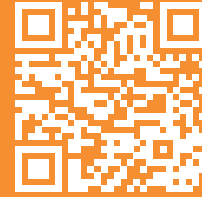

INDIVIDUAL ORDER  
NUMBERS CAN BE FOUND AT  
[RELIANCEORTHODONTICS.COM](http://RELIANCEORTHODONTICS.COM)

### ✓ NITI & SS Stainless Steel 25/PKG TM Beta Titanium 10/Pkg

✓ = AVAILABLE IN ALL TNT ARCHFORMS

ALLOYS	NICKEL TITANIUM		TM BETA TITANIUM	SS STAINLESS STEEL
	HEATERS THERMAL NITI	SE SUPERELASTIC NITI		
.012	✓	✓	-	✓
.013	✓	✓	-	-
.014	✓	✓	-	✓
.016	✓	✓	-	✓
.018	✓	✓	-	✓
.020	✓	✓	-	✓
.014 x .025	✓	✓	-	-
.016 x .016	✓	✓	-	✓
.016 x .022	✓	✓	✓	✓
.016 x .025	✓	✓	-	✓
.017 x .025	✓	✓	✓	✓
.018 x .018	✓	✓	-	✓
.018 x .025	✓	✓	✓	✓
.019 x .025	✓	✓	✓	✓
.020 x .020	✓	✓	-	✓
.021 x .025	✓	✓	✓	✓

### ✓ Specialty NITI Archwires



10/PKG

#### TNT REVERSE CURVE NITI

R-FORM UP & LOW D-FORM UP

(Available only in SE NITI sizes)

.014	.014 x .025	.016 x .016
.016	.016 x .022	.016 x .025
.018	.017 x .025	.019 x .025
.020	.021 x .025	




25/PKG

#### TNT PRESTOPPED HEATERS, SE OR SS

ALL TNT ARCHFORMS

(Not available in .021 x .025 wires)

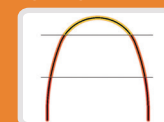


10/PKG

#### TNT PROGRESSIVE FORCE NITI

R-FORM UP & LOW

(Available in all dimensions)

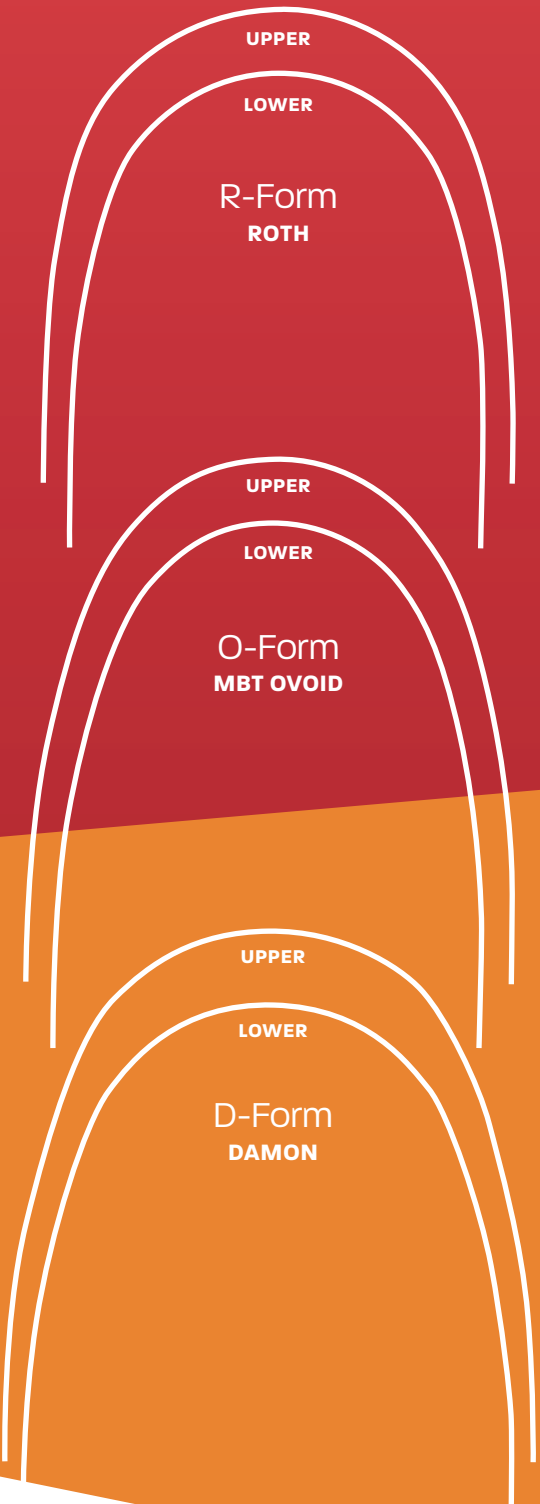

# Relastics



## A New Standard of Consistency in Intraoral Forces

Reliance, the name you trust in bonding, has researched materials and manufacturing processes to develop a precision line of patient elastics. The result is unprecedented consistency and accuracy in forces as prescribed.

### Relastics Only From Reliance

- 5 sizes and 4 forces provide the treatment choices for any applied mechanics
- Manufactured from surgical grade latex
- Thickness and diameters guaranteed accurate
- Precision manufacturing techniques and in-process testing assure force accuracy
- Force values determined at 3 times extension from passive diameter
- Patient pack ID marks are familiar and concise
  - Letters designate size
  - Colors designate force
- Light retardant packaging to eliminate UV degradation
- Resealable and tamper-proof patient packs
- 3-year shelf life



- ✓ Reliance Orthodontic Products Over 40 Years of Serving the Orthodontic Profession.



Letter & Color IDs



Resealable & Tamper Proof

Letter Designates Size

Color Designates Force

Light Resistant Plastic

Relastics Size Guide

50 Packs of 100 Elastics Per Pouch



Force Values & Order Numbers

Actual Size	Dimensional Size	Size ID	Light 2.5 oz Red	Medium 3.5 oz Green	Heavy 4.5 oz Blue	X Heavy 6.5 oz Black
	1/8" 3.2mm	A	Red A #EL18L	Green A #EL18M	Blue A #EL18H	Black A #EL18X
	3/16" 4.8 mm	B	Red B #EL36L	Green B #EL36M	Blue B #EL36H	Black B #EL36X
	1/4" 6.4 mm	C	Red C #EL14L	Green C #EL14M	Blue C #EL14H	Black C #EL14X
	5/16" 7.9 mm	D	Red D #EL56L	Green D #EL56M	Blue D #EL56H	Black D #EL56X
	3/8" 9.5 mm	E	Red E #EL38L	Green E #EL38M	Blue E #EL38H	Black E #EL38X

# Custom Base Indirect Bonding

## MODEL ADHESIVES



- Tacky paste viscosity
- Unlimited working time
- Heat curable
- Highly filled composite



- Viscous to resist drift on resin and stone models
- Extended working time to allow bracket manipulation
- Single-piece flash clean-up

THERMA-CURE™	
5g Therma-Cure™ Syringe	#TC
GO TO®	
4g Paste Syringe	#GTP
4g Paste In Tips (10 Tips)	#GTT

## TRANSFER ADHESIVES



- Needle Tip Dispenser
  - Using the tips provided, these materials are easily injected onto custom pads
- Low modulus of elasticity to resist breakage
- Highly filled to bridge gaps between custom pads and tooth surface



- 90-Second working time
- Excellent molecular conversion rate
- Self-curing

FLOWTAIN™	
1.5g Flowtain™ with Tips	#FT
19 gauge Tips (10)	#19TIPS10
MAXIMUM CURE® SEALANT	
7g Maximum Cure® Sealant A with Fluoride	#MC7AF
7g Maximum Cure® Sealant B	#MC7B
7g Maximum Cure® Filled Sealant A with Fluoride	#MCF7AF
7g Maximum Cure® Filled Sealant B	#MCFSB
20g Maximum Cure® Sealant A with Fluoride	#MC20AF
20g Maximum Cure® Sealant B	#MC20B

# Clean Base IDB Adhesives

**#1** OPTION FOR DIGITALLY GUIDED PROCESS

**L.C.R.™**  
HYBRID FLOWABLE IDB ADHESIVE



L.C.R.™	
5g L.C.R.™ Paste in Push Syringe	#LCR
1.5g L.C.R.™ Paste in Push Syringe w/ Tips	#LCR15
5g L.C.R.™ Paste in Tips (11)	##LCRT5

- Hybrid flowable consistency allows for instant bracket mesh penetration
- Resin matrix holding capacity allows for preloading up to 7 days in advance
- Micro+ nano fill for flowable consistency with durability of highly filled composite

**GO TO®**  
LIGHT CURE BRACKET ADHESIVE



GO TO®	
4g Paste Syringe	#GTP
4g Paste In Tips (10 Tips)	#GTT

- Perfect viscosity to allow bracket mesh penetration in jigs
- Highly filled to provide ultimate strength
- Pre-loadable – allowing even shorter chair time prep
- Fluoresces for easy removal

**EXCEL®**  
CHEMICAL CURE SPECIALTY ADHESIVE



EXCEL®	
Excel® Regular Paste A in 6g Screw Syringe with Fluoride	#1ERSAF
Excel® Regular Paste B in 6g Screw Syringe	#1ERSB
Excel Regular Paste A in 12g Jar with Fluoride	#ERJAH
Excel Regular Paste B in 12g Jar with Fluoride	#ERJB

- Unsurpassed bond strength
- Thin, creamy consistency
- Facilitates clean base penetration
- Three minute working time
- Contains fluoride

\* Excel® Paste available w/o Fluoride on request

# Accessories

## DISPOSABLE BENDA BRUSHES



DISPOSABLE BENDA BRUSHES

### DISPOSABLE BENDA BRUSHES

4" Benda Brushes (250)	#4BB250
------------------------	---------

- Sealant and enhancer applicators

## NYLON MIXING SPATULAS



NYLON MIXING SPATULAS

### NYLON MIXING SPATULAS

Spatulas (50)	#S
---------------	----

## MICRO BRUSHES



MICRO BRUSHES

### MICRO BRUSHES

Micro Brushes (100)	#MB
Super Fine Micro Brushes (100)	#MBSF

## SPONGE PELLETS



SPONGE PELLETS

## BRUSH HANDLE AND TIPS



BRUSH HANDLE AND TIPS

### SPONGE PELLETS

Sponge Pellets (500)	#SP
----------------------	-----

### BRUSH HANDLE AND TIPS

Brush Handle	#BHA
Brush Tips (250)	#BT250

- Etchant applicators

 MIXING PADS

 MIXING WELL




MIXING PADS

Disposable Mixing Pads (1) #1PADS

MIXING WELL

Mixing Well #1MXW

 S.E.P. MIXING WELL

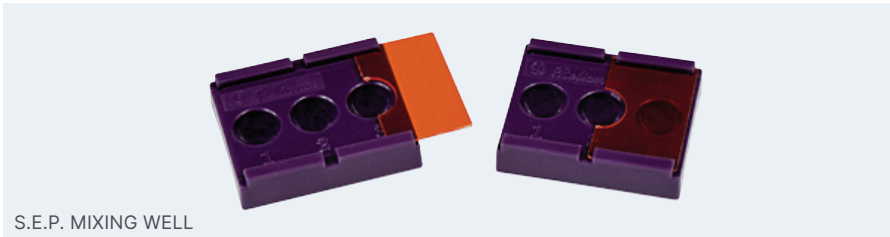

S.E.P. MIXING WELL

S.E.P. Mixing Well with Light Impervious Shield #1SEPW

- Ambient light blocking

 UVEX EYEWEAR



UVEX EYEWEAR

Uvex Ultra Spec® 2000 #SGU

- Fits over prescription eyeglasses

 SKYPER® EYEWEAR



UVEX SKYPER®

UVEX skyper® #SG

- Low profile

# Accessories

## MARK II CR SYRINGE



MARK II CR SYRINGE


MARK II CR SYRINGE

Mark II CR Syringe

#CRII

- Used with pre-loaded adhesive tips

## DIAMOND BUR 129

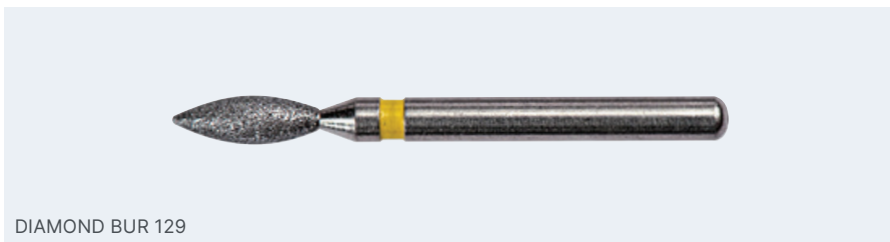

DIAMOND BUR 129


DIAMOND BUR 129

Diamond Bur 129

#DB

- Prep bur
- Roughens surface
- \* Other burs appear on page 23


## ACCESS™ RETRACTOR

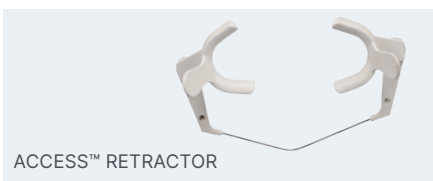

ACCESS™ RETRACTOR

- Sterilize with dry heat or autoclave

## TECHNIQUE CARDS

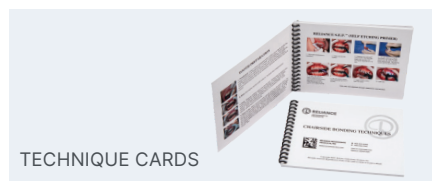

TECHNIQUE CARDS

- Chairside reference for all bonding scenarios
- This item is non-returnable

ACCESS™ TRACTOR


Access™ Adult Cheek Retractor #CH1

Access™ Pado Cheek Retractor #CH2

TECHNIQUE CARDS

Chairside Bonding Technique Cards #TECH

## KOBAYASHI TIES

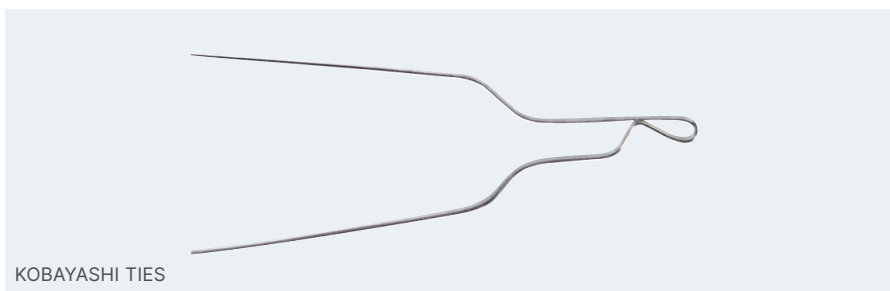

KOBAYASHI TIES

KOBAYASHI TIES

Tie Hooks - .010 100/Pkg #K10

Tie Hooks - .012 100/Pkg #K12

Tie Hooks - .014 100/Pkg #K14

Pre-Twisted - .012 100/Pkg #SK12

Pre-Twisted - .014 100/Pkg #SK14

# Terms & Conditions

---

## ORDERING

To order Reliance Orthodontic Products, simply contact our home office Toll-Free (in U.S. and Canada) at 1-800-323-4348 or 1-630-773-4009. You can also order through our secure website at [www.RelianceOrthodontics.com](http://www.RelianceOrthodontics.com). **There is a minimum of \$25.00 on all orders.**

## SHIPPING

All prices are F.O.B. Itasca and are shipped via FedEx. All shipping is billed with the order. Alternate shipping services are available. Please call for more information.

## EXPORT

Export orders must be paid for in advance by proforma invoice.

## TERMS

Net 30 days. Illinois residents please add 1.75% sales tax. An invoice is enclosed with each shipment and payment should be made from this invoice. Any outstanding balance will be subject to a finance charge of 1.5% per month. All payments are due in U.S. funds.

## EXCHANGES OR REFUNDS

Your satisfaction is Guaranteed! Reliance Orthodontic Products, Inc. guarantees that our products are of the highest quality and workmanship. If at any time prior to its expiration date your adhesive system fails to perform in a manner not consistent with the product, we will replace the remaining product with new material free of charge or issue a credit for the remaining product. Please check your order for accuracy. Any discrepancies must be reported to Reliance within 7 days. No refunds, credits, or replacement will be issued after 30 days unless the product is deemed defective by Reliance Orthodontic Products, Inc. Products that are being returned to Reliance Orthodontic Products, Inc. must have a prior Return Authorization Number (RMA). All products being returned for credit or replacement must be sent postage paid with a copy of your invoice and your RMA number clearly visible on the package, or the package will be refused by Reliance Orthodontic Products, Inc. A restocking fee may be applicable.

## CONTACT

1540 West Thorndale Ave., PO Box 678,  
Itasca, IL 60143 USA

1-800-323-4348

Phone 1-630-773-4009

Fax 1-630-250-7704

Please email any questions or comments to [pgange@RelianceOrthodontics.com](mailto:pgange@RelianceOrthodontics.com) or [paulgangejr@RelianceOrthodontics.com](mailto:paulgangejr@RelianceOrthodontics.com)



RELIANCEORTHODONTICS.COM

© 2022 Reliance Orthodontic Products




US FDA  
**REGISTERED**



HEALTH CANADA  
**REGISTERED**

Certificates Available Upon Request

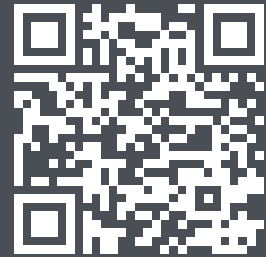

MADE IN THE U.S.A.